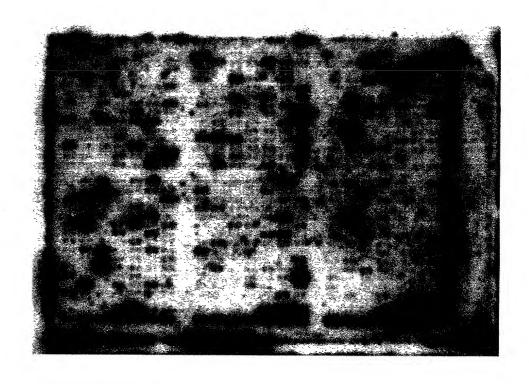
Tissue or cells • Isolate total RNA • Treat with DNase I Immit Poly A⁺ RNAs anneal with Biotinylated Oligo(dT) Biotinylated Oligo(dT) Streptavidin Magnetic Beads · Biotin binds to Streptavidin Magnetic Beads • Wash away poly ATRNAs and other impurities Add gene-specific primer mix & other reagents Atlas™ cDNA probe synthesis provided with your Atlas™ Array and labeling procedure

Elute probe from beads, purify, and hybridize to Atlas™ Array



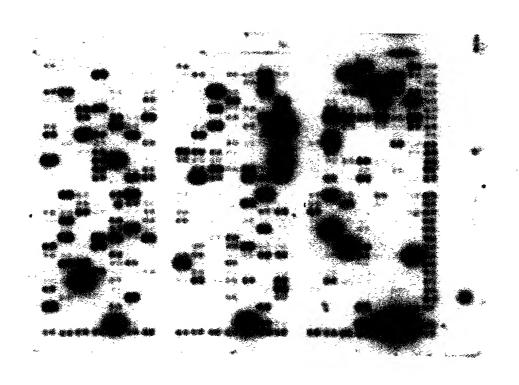


FIGURE 2